

260  
A.O  
W2  
MA  
S9  
ITEM NO. 22  
FILE NO. XXVIII-33

QV  
354  
P337r  
1947

COPY NO. 181

1947  
NO. 28-33  
1947  
UNCLASSIFIED

CONFIDENTIAL

CLASSIFICATION ORIGIN

UNCLASSIFIED



Sec I Vol I, J10A, Acc List

26 June 1947

*Frank B Rogers*

# EXPERIMENTAL PRODUCTION OF PENICILLIN AT THE ST. JACOB'S HOSPITAL, LEIPZIG

Exec	Prof	.....
JSS	Med	.....
SBJ	NP	.....
ESR <i>Col. Long</i>	Surg	.....
Trop	Recon	.....
VD	Vet	.....
Lab	Oper	.....
San	Train	.....
SE	Hosp	.....
Civ-Nutr	Supp	.....
.....	Control	.....
.....	Med Stat	.....
.....	Hist	.....
.....	Spec Dev	.....

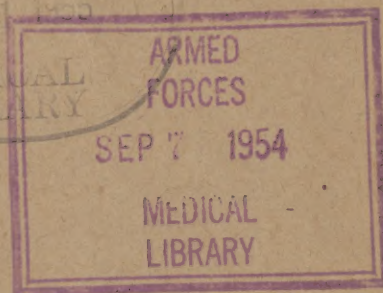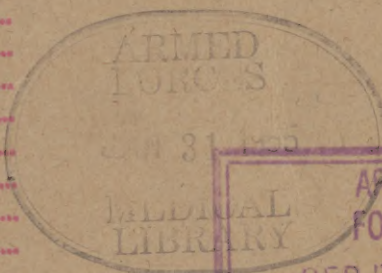

*Army Medical School*  
Initial and Return to  
Medical Intelligence  
NOT  
to next office checked

CONFIDENTIAL

UNCLASSIFIED

COMBINED INTELLIGENCE OBJECTIVES  
SUB-COMMITTEE





CONFIDENTIAL

Report on

EXPERIMENTAL PRODUCTION OF PENICILLIN AT THE ST. JACOB'S  
HOSPITAL, LEIPZIG

Reported by:

Mr. P.L. PAVCEK, T.I.I.C

20097

CIOS Target No. 22/1289  
Food and Agriculture

COMBINED INTELLIGENCE OBJECTIVES SUB-COMMITTEE  
G-2 Division, SHAEF (Rear) APO 413

CONFIDENTIAL



## ST. JACOB'S HOSPITAL.

### Department of Dermatology.

-----

Dr. Joseph Kimmig started research on penicillin under Prof. Joseph von Kennel in Kiel in 1942. In 1943 the firm of Schering A.G., supporting the work, transferred the project to St. Jacob's Hospital where it has been in progress ever since. The Hospital is a very new building and has beautiful laboratories and equipment.

Penicillin was being grown on a rather small scale in one of the laboratories in bottles containing 300 ccs. of Chopek-Dox medium. No other adjuncts were added. The mold, of which a culture was obtained, was isolated by Kimmig from the flower of Blue Berries and was identified as P. Chryseum flavum. As high as 40 units of Penicillin per cc. were produced in a ten day incubation at 26°C. but for continuous operation the figure was 10-15 units/cc.

A special procedure was adopted to build up the inoculum. The spore suspension grown on a standard nutrient Kuczynski agar (sample obtained) was diluted and plated out until a culture originating from 1 spore was obtained. 12-15 plates of such spore cultures were investigated separately for ability to produce Penicillin, and in this manner the most potent culture was selected.

For testing the unitage in the final solution, 1 cc. of the test solution was inoculated with S. aureus incubated at 37° for 24 hours. Dilutions were made according to the attached protocol (Table 1).

The total weekly capacity of the laboratory set-up was 200,000 units. Some experiments had been conducted with submerged cultures but maximum unitage obtained was 1/4 Unit/cc. Kimmig felt that this culture was not suited to this type of culturing.

Other workers on Penicillin in Germany included: Bernhauer in Prague, R. Kuhn (Heidelberg) and the firm Schatt & Genossen (Jena).

Dr. Kimmig demonstrated a new observation on inhibition of oxygen uptake by S. aureus in the presence of Penicillin or malonic acid. On a Ringer phosphate medium 1:10,000 malonic acid produced the same inhibition

as did Penicillin. According to Kimmig the effect of Penicillin may be due to a disturbance in the succinic acid cycle.

Dr. Hermann Fox also associated with Schering A.G. had done research on sulfonamides.

A sulfonamide like compound had been produced by I.G. Farben (Leverkusen) and had very powerful bactericidal effects. The compound, tetra-bromo-benzil, had not been clinically tested as yet.

Several reprints from Kimmig were obtained. One on the effect of para amino benzo acid in preventing sunburn was especially interesting.



# Testschema für Kulturflüssigkeit

Datum: \_\_\_\_\_

Fraktion: \_\_\_\_\_

Verdünnungsstufe	ccm	Wachstum	PE
VSt D I	0.1 + 5 Standard (1 × konz.)		500
2 ccm D + 8 ccm	0.15 +		330
Puffer (1 : 500)	0.2 +		250
	0.25 +		200
	0.3 +		170
	0.35 +		140
	0.4 +		125
	0.45 +		110
	0.5 +		100
VSt D	0.1 +		100
1 ccm C + 9 ccm	0.2 +		50
Puffer (1 : 100)	0.25 +		40
	0.3 +		33
	0.4 +		25
	0.5 +		20
VSt C I	0.2 +		25
2 ccm C + 8 ccm	0.25 +		20
Puffer (1 : 50)	0.3 +		17
	0.4 +		13
	0.5 +		10
VSt C	0.1 +		10
1 ccm Kfl + 9 ccm	0.13 +		8
Puffer (1 : 10)	0.2 +		5
	0.3 +		3.5
	0.5 +		2
VSt B	0.1 +		1
Unverd. Kfl	0.2 +		0.5
	0.3 +		0.35
	0.5 +		0.2

## Kontrolle ohne Zusatz

### VSt A



ccm Kfl	ccm Puffer	ccm 2 × konz. Std.	Wachstum	PE
0.5	2	2.5		0.2
1.0	1.5	2.5		0.1
1.5 ✓	1.0	2.5		0.07
2 ✓	0.5	2.5		0.05
2.5 ✓	0	2.5		0.04





Frühwasser von Ges. 9 25°C  
mit 17 cm Wasserdruck

Frühwasser von Ges. 9 mit Kesselwasser-  
druck über 12,8 bar Vorwärmer & Umwälzpumpe 25°C

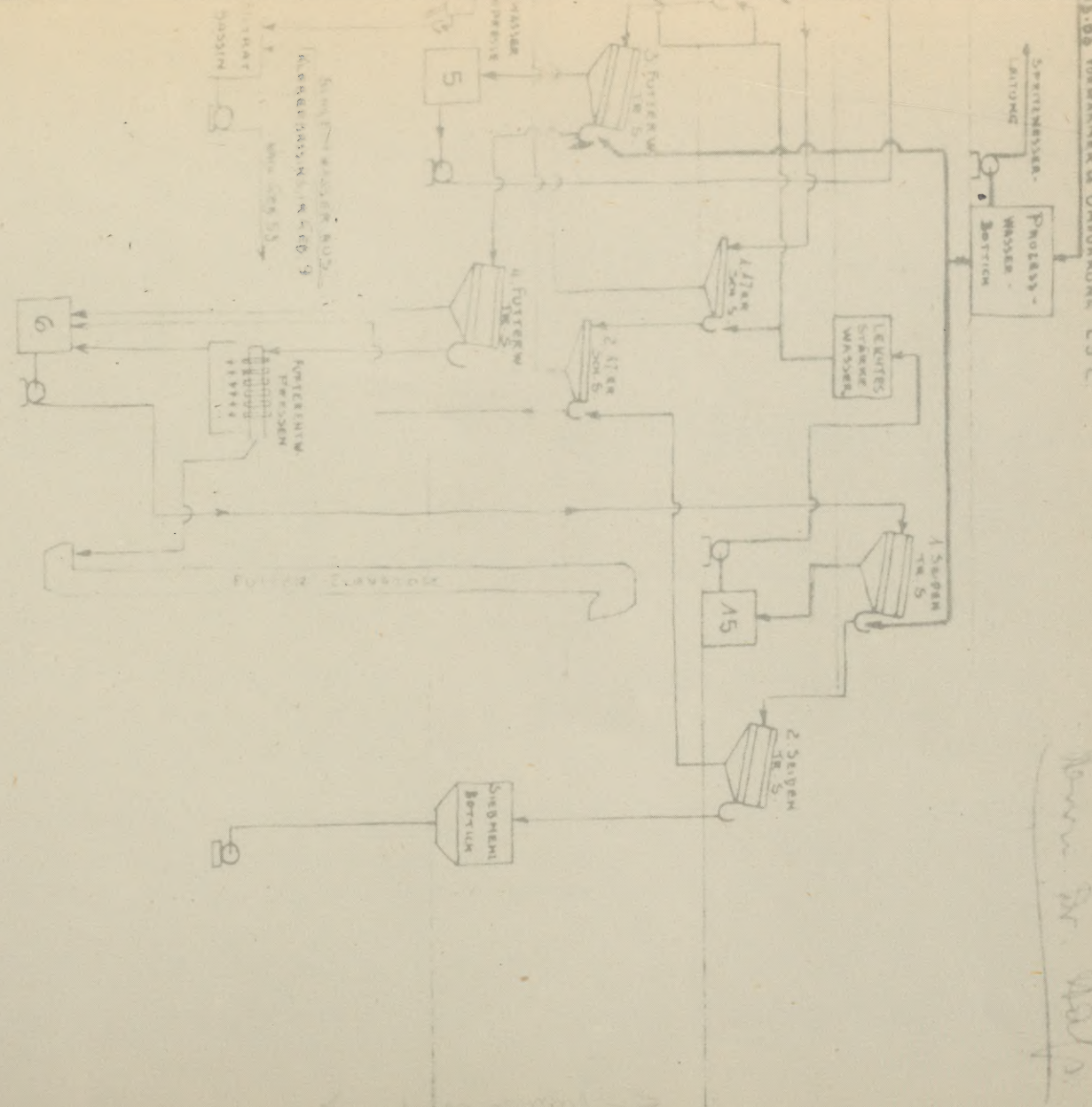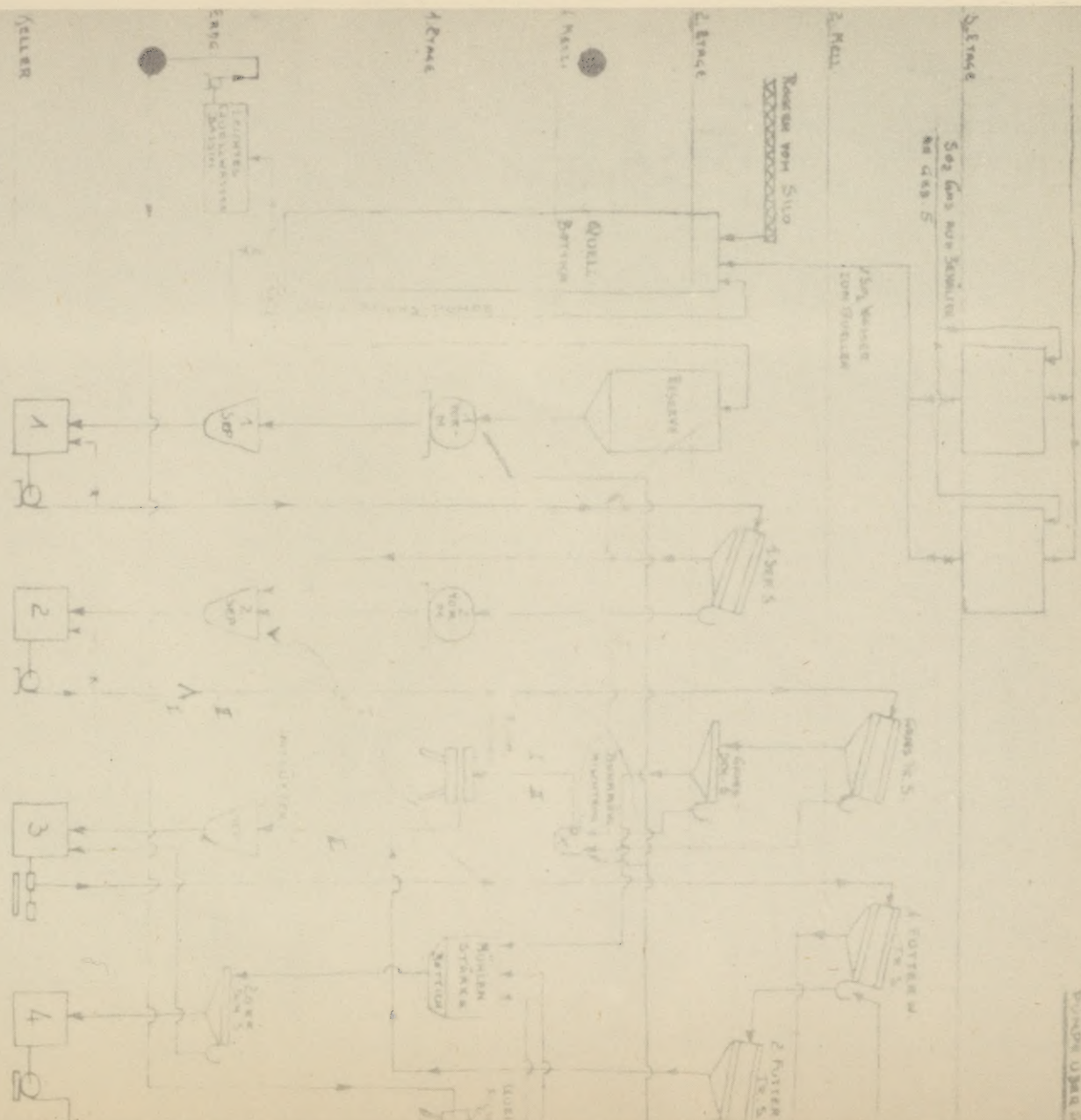

# ARBEITS-SCHEMA

FÜR DIE ROGGENVERARBEITUNG VOM

8 MÄRZ BIS

1945

BARBY/ELBE, 8. MÄRZ 45 MÜ.

FIGURE 1. FLOW DIAGRAM OF PLANT FOR STRACH  
MANUFACTURE FROM RYE

Wasser 17 cm  
alle anderen  
12,8 bar  
17 cm  
12,8 bar

Wasser 17 cm











NATIONAL LIBRARY OF MEDICINE



NLM 03617446 3

20097